


## DL700 Digital Holter Monitor

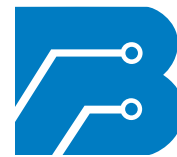

**Braemar, Inc.**

### Digital Holter Monitor

- Flash Memory
  - Nonvolatile
  - No data compression
- Precise Digital Recording
- PCMCIA Compatibility – (Media Plugs into PC)
- 24-Hour Battery Operation
- 3-Channel Recording
  - 7 wire patient interface
- Proven, Reliable and Simple Operation
- Patient Event Button
- LCD Clock
- AA Battery Operation

*The Braemar DL700 Digital Holter is a highly-reliable monitor designed for continuous 24-hour recording. Its long history of performance ensures superior reliability and data integrity while maximizing patient convenience and safety.*

*The flash memory card with patient data can be easily removed by “ejector button,” freeing the monitor to be used by another patient while the data is loaded via an industry-standard PCMCIA slot.*

*The DL700 is designed for custom configuration to provide unique functionality for volume OEM applications. A recessed area on the top cover provides space for your company identification.*

# DL700

## Digital Holter Monitor

### Product Specifications

#### FUNCTIONAL

Channels ..... 3

#### MEMORY

Recording ..... Full Disclosure  
24 Hours

Flash Memory ..... PCMCIA Flash Card, Nonvolatile

#### PHYSICAL

Dimensions ..... 6" x 3.5" x .95"  
(152mm x 89mm x 24mm)

Weight ..... 7.6 oz (215 grams)

Enclosure ..... Molded Plastic / Water Resistant

#### ELECTRICAL

Connector ..... PCMCIA-ATA I/O Nonvolatile

Resolution ..... 8 Bit, 128/133 Samples Per Second

Test Jacks ..... 3 Channel, Color Coded

Frequency Response ..... 0.05 - 60Hz

Patient Connection ..... 7 Pin Patient Wire, Direct Connect,  
Color Coded, Snap On

#### BATTERY

Type ..... (2) AA Alkaline

**WARRANTY** ..... 1 Year

#### OEM OPTIONS

- Cover Colors: White Standard, Others Optional
- Patient Cable
- Custom Configurable
- Flash Card Data Capacity
- PC Data Transfer Utility Software



**Braemar, Inc.**  
1285 Corporate Center Drive  
Suite 150  
Eagan, MN 55121 USA

[www.braemarinc.com](http://www.braemarinc.com)

[sales@braemarinc.com](mailto:sales@braemarinc.com)

800.328.2719

651.286.8620

651.286.8630 fax